



# Increase Monoclonal Antibody Purity with TOYOPEARL® Sulfate-650F C-IEX Resin for Capture and Removal of mAb Aggregates

## INTRODUCTION

Ion exchange chromatography is often used as an intermediate purification step in monoclonal antibody (mAb) purification for the removal of protein aggregates, host cell proteins (HCP), and leached protein A ligand. Industry trends are focusing on the development of continuous downstream processing. Typically scientists in biopharmaceutical settings use cation exchange (CEX) and anion exchange (AEX) chromatography steps in series to further polish a purified mAb after the protein A capturing step. In this study, we focus on the development of a post protein A CEX step for the removal of aggregates, host cell proteins and leached protein A to improve the purity of the mAb eluate in a single polishing step.

A strong cation exchange resin, TOYOPEARL Sulfate-650F, is used in this study. It is a novel resin with the following benefits: strong capture of mAb aggregates, high salt-tolerance, wide working pH range, and high dynamic binding capacity. The protocol used with this resin has a minimal pH adjustment of the purified mAb eluate after the protein A step.

## MATERIALS AND METHODS

TOYOPEARL Sulfate-650F resin, 45  $\mu\text{m}$ , 100 nm and TOYOPEARL AF-rProtein A HC-650F, 45  $\mu\text{m}$ , 100 nm resins were used in this study. The resins were packed into Omnifit® Benchmark columns, (25 mm ID  $\times$  5 cm or 6.6 mm ID  $\times$  10 cm). A TSKgel G3000SW<sub>XL</sub> SEC column was used for analyzing collected IgG<sub>1</sub> fractions

### Purify IgG<sub>1</sub> using TOYOPEARL AF-rProtein A HC-650F

As shown in the steps below, CHO clarified cell supernatant (CCS) containing IgG<sub>1</sub> was passed through a TOYOPEARL AF-rProtein A HC-650F column (25 mm ID  $\times$  15 cm) to purify the IgG<sub>1</sub>. The elution peak was collected starting and ending at 100 mAU.

The IgG<sub>1</sub> eluate from the TOYOPEARL AF-rProtein A HC-650F column was adjusted to pH 5.0 with 1 mol/L Tris base and quantified by UV absorbance at 280 nm.

- 1 Equilibrate (6 CV, 225 cm/h): 20 mmol/L Tris-acetate, 150 mmol/L NaCl, pH 7.4
- 2 Load (48 mg/mL-resin, 225 cm/h): TBL-mAb-01 CCS
- 3 Wash (10 CV, 225 cm/h): 20 mmol/L Tris-acetate, 150 mmol/L NaCl, pH 7.4
- 4 Elute (5 CV, 225 cm/h): 50 mmol/L acetic acid, pH 3.0
- 5 Sanitize (4 CV, 225 cm/h): 0.1 mol/L NaOH
- 6 Regenerate (4 CV, 225 cm/h): 20 mmol/L Tris-acetate, 150 mmol/L NaCl, pH 7.4

### Static Binding Capacity Screening

To obtain static binding capacity (SBC) for the mAb on TOYOPEARL Sulfate-650F, the collected purified IgG<sub>1</sub> eluate was dialyzed and adjusted to various pH and sodium chloride conditions and a final total IgG<sub>1</sub> concentration of 10 mg/mL. The IgG<sub>1</sub> was bound to TOYOPEARL Sulfate-650F resin in Resin Seeker plates (20  $\mu\text{L}$ /well) as described below. Adjusted protein A IgG<sub>1</sub> eluate was bound to the resin by batch adsorption for 1 hour at ambient temperature. Following incubation, resin was removed from each well by vacuum filtration, and 75  $\mu\text{L}$  samples of each well were read for UV absorbance to determine concentration of unbound protein. Static binding capacity was determined. SBC data was analyzed with SAS JMP 12 software.

- 1 Equilibrate (3  $\times$  200  $\mu\text{L}$ /well): 50 mmol/L Tris-acetate\*, various NaCl conc. and pH
- 2 Bind (200  $\mu\text{L}$ /well, 1 h, RT) : TBL-mAb-01, 10 g/L, various NaCl conc. and pH

\*Note: for pH 4.0 – 5.6, 50 mmol/L acetic acid was titrated with Tris base. For pH 7.2 – 8.4, 50 mmol/L Tris base was titrated with acetic acid.

### DYNAMIC BINDING CAPACITY OPTIMIZATION

Dynamic binding capacity (DBC) at 10% breakthrough was determined for TOYOPEARL Sulfate-650F. Dialyzed protein A eluate was adjusted to various pH and sodium chloride concentrations and a final IgG<sub>1</sub> concentration of 5 mg/mL. DBC determination was done as shown below. Chromatograph was primed with protein solution to determine UV absorbance (280 nm) at 100% breakthrough.

Protein was bound to column until UV absorbance at 10% breakthrough was reached, and DBC was determined based on the volume of protein solution loaded. DBC data was analyzed with SAS JMP 12 statistical software.

- 1 Equilibrate (10 CV, 180 cm/h): 50 mmol/L acetate-Tris, various NaCl conc. and pH
- 2 Load (45 cm/h): TBL-mAb-01, 5 g/L to ca. 136 mAU
- 3 Wash (5 CV, 45 cm/h): equilibration buffer
- 4 Elute (5 CV, 45 cm/h): equilibration buffer + 1 mol/L NaCl
- 5 Sanitize (5 CV, 60 cm/h): 0.5 mol/L NaOH
- 6 Regenerate (5 CV, 60 cm/h): water

### ELUTION OPTIMIZATION

Collected IgG<sub>1</sub> eluate was adjusted to pH 5.2 with 1 mol/L acetic acid and/or 1 mol/L Tris base and 12.1 mS/cm conductivity with 4 mol/L NaCl and/or water. Sample was loaded onto a 6.6 mm ID × 3.0 cm column of TOYOPEARL Sulfate-650F as shown below. Conductivity at elution peak was determined.

- 1 Equilibrate (10 CV, 180 cm/h): 50 mmol/L acetate-Tris, 100 mmol/L NaCl, pH 5.2
- 2 Load (45 cm/h): TBL-mAb-01, ca. 10 g/L
- 3 Wash (5 CV, 45 cm/h): equilibration buffer
- 4 Elute (20 CV, 45 cm/h): 50 mmol/L acetate-Tris, 100 – 500 mmol/L NaCl, pH 5.2
- 5 Sanitize (5 CV, 60 cm/h): 0.5 mol/L NaOH
- 7 Regenerate: (5 CV, 60 cm/h): water

Separation was repeated with a 98 mg/mL-resin load, and step gradient elution at 260, 290, or 320 mmol/L NaCl for 10 CV. Throughout step gradient elution, 1-CV fractions were collected and analyzed for protein concentration, aggregate content (SEC chromatography with TSKgel G3000SW<sub>XL</sub>), CHO-HCP (ELISA), and protein A content (ELISA).

- 1 Equilibrate (10 CV, 180 cm/hr): 50 mmol/L acetate-Tris, 100 mmol/L NaCl, pH 5.2
- 2 Load (45 cm/h): TBL-mAb-01, ca. 10 g/L
- 3 Wash (5 CV, 45 cm/h): equilibration buffer
- 4 Elute (10 CV, 45 cm/h): equilibration buffer + 260 – 320 mmol/L NaCl
- 5 Strip (5 CV, 45 cm/h): 50 mmol/L acetate-Tris, 1.0 mol/L NaCl, pH 5.2
- 6 Sanitize (5 CV, 60 cm/h): 0.5 mol/L NaOH
- 7 Regenerate: water (5 CV, 60 cm/h)

### RESULTS AND DISCUSSIONS

#### Purification of IgG<sub>1</sub> with TOYOPEARL AF-rProtein A HC

The crude sample containing IgG<sub>1</sub> was passed through the protein A column and fractions of IgG<sub>1</sub> were collected for further work. Figure 1 demonstrates that the IgG<sub>1</sub> was purified by protein A chromatography. The eluate peak was collected and further analyzed by size exclusion chromatography using a TSKgel G3000SW<sub>XL</sub> SEC column for monomer and aggregate yield, host cell protein (HCP) content

#### PURIFICATION OF IgG<sub>1</sub> FROM CHO SUPERNATANT

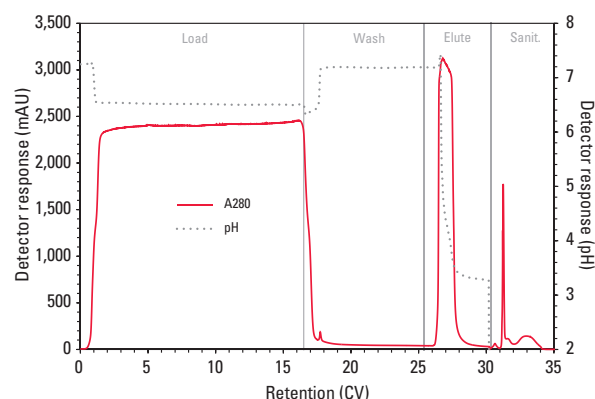


Figure 1

Resin: TOYOPEARL AF-rProtein A HC-650F  
 Column: 25 mm ID × 15 cm (74 mL)  
 Mobile phase: A: 20 mmol/L Tris-acetate, 150 mmol/L NaCl, pH 7.4  
 B: 50 mmol/L acetic acid  
 C: 0.1 mol/L NaOH  
 Flow rate: 225 cm/h (4 min residence time)  
 Detection: UV @ 280 nm (mAU), pH  
 Temperature: ambient  
 Injection vol.: 1200 mL (48 mg/mL-resin load ratio)  
 Sample: TBL-mAb-01 CSS @ 2.95 g/L and protein A ligand leaching (see Table 1).

#### ANALYSIS DATA FOR THE COLLECTED IgG<sub>1</sub> ELUATE PEAK

Protein A Eluate Analysis	
Yield (total IgG)	99%
Aggregate	4.4% (0.5% HMW, 3.9% dimer)
HCP	1260 ppm
Protein A	1.2 ppm

Table 1

#### Static binding capacity (SBC) screening for TOYOPEARL Sulfate-650F

To optimize the binding capacity of TOYOPEARL Sulfate-650F resin, the SBC screening was carried out to find out the maximum amount of protein bound to a chromatography medium at given solvent and protein concentration conditions. Fig. 2 shows a binding map that was created for TOYOPEARL Sulfate-650F resin at pH values from 4.0 – 5.6 or 7.2 – 8.4, and 0 – 500 mmol/L NaCl. Maximum protein binding was noted between pH 4.8 and 5.6, from ca. 100 mmol/L – 200 mmol/L NaCl. At lower pH values, NaCl is necessary for protein binding. At higher pH values, significant binding is noted with little to no NaCl.

#### Dynamic Binding Capacity (DBC) Optimization

DBC was optimized by DoE with a three-level, full-factorial method at pH 4.8 – 5.6 and 100 – 200 mmol/L NaCl (results as shown in Figure 3 with data points consolidated in Table 2). A maximum DBC of >120 mg/mL-resin was noted between pH 4.8, 150 mmol/L NaCl, and pH 5.2, 100 mmol/L NaCl. Conditions of pH 5.2, 12.1 mS/cm were used for the elution optimization experiments for maximum binding.

SBC SCREENING FOR TOYOPEARL SULFATE-650F USING THE COLLECTED IgG<sub>1</sub> ELUATE

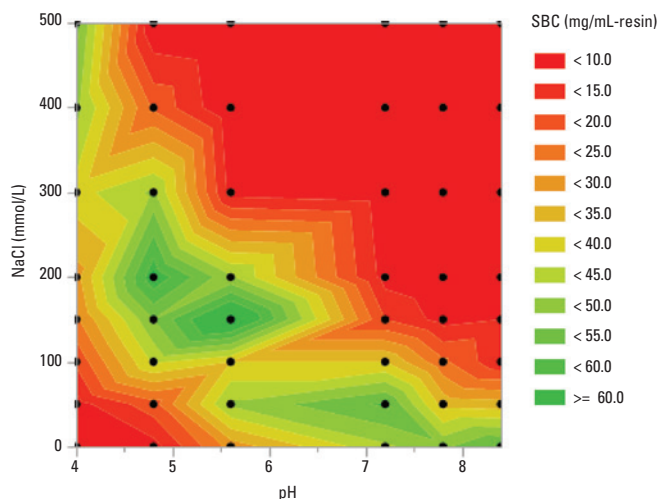


Figure 2

DYNAMIC BINDING CAPACITY OPTIMIZATION FOR TOYOPEARL SULFATE-650F USING THE COLLECTED IgG<sub>1</sub>

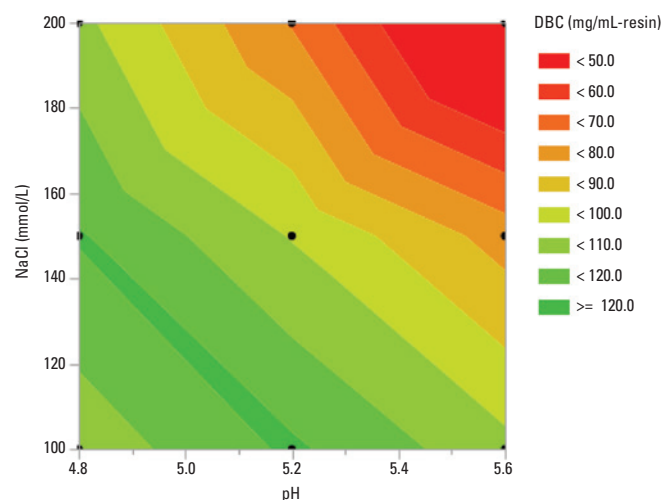


Figure 3

DYNAMIC BINDING CAPACITY DATA AT VARIOUS CONDITIONS (RUN ORDER OF CONDITIONS RANDOMIZED)

Exp.	Run	Load pH	Load NaCl (mmol/L)	DBC (mg/mL)
1	5	4.8	100	104
2	1	5.6	100	103
3	2	4.8	200	103
4	4	5.6	200	23
5	8	4.8	150	121
6	9	5.6	150	76
7	7	5.2	100	122
8	6	5.2	200	69
9	3	5.2	150	99
10	10	5.2	150	102

Table 2

Optimize conditions for separation and elution of monomer and high molecular weight peaks

To optimize elution conditions, a gradient elution was performed at pH 5.2. A peak conductivity of 30.1 mS/cm was noted (ca. 288 mmol/L NaCl). Experiment was repeated as a step gradient at 260, 290, or 320 mmol/L NaCl (see Fig 4). Due to peak tailing during elution, 1-CV fractions were collected throughout elution. Fractions were analyzed for IgG<sub>1</sub> concentration, aggregate, HCP and protein A, and results were analyzed to determine optimum NaCl concentration and peak volume.

Peaks were analyzed for recovery, aggregates, HCP and protein A content (Figure 5). Data analysis suggests the optimum aggregate and HCP removal are obtained at 260 mmol/L NaCl in elution buffer and maximum (9 CV) elution volume. Protein A ligand content at these conditions (40 ppb) is significantly lower than that found in the load material (1200 ppb). Data is consolidated in Table 3.

PROFILING OF THE COLLECTED IgG<sub>1</sub> ELUATE PEAK SEPARATED BY TOYOPEARL SULFATE-650F AT VARIOUS CONDITIONS

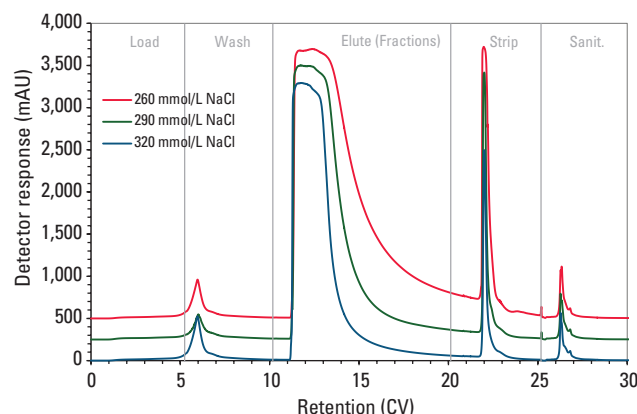


Figure 4

Resin: TOYOPEARL Sulfate-650F Column; 6.6 mm ID x 3.0 cm (1.0 mL)  
 Mobile phase: A: 50 mmol/L acetate-Tris, 100 mmol/L NaCl, pH 5.2  
 B: 50 mmol/L acetate-Tris, pH 5.2, NaCl as indicated  
 C: 50 mmol/L acetate-Tris, 1.0 mol/L NaCl, pH 5.2  
 D: 0.5 mol/L NaOH  
 Flow rate: 45 cm/h (4 min residence time)  
 Detection: UV @ 280 nm (mAU)  
 Temperature: ambient  
 Injection vol.: 5.3 mL (97 mg/mL-resin load ratio)  
 Sample: TBL-mAb-01, 19.1 mg/mL

Exp.	Run	Elution NaCl (mmol/L)	Pool volume (CV)	Recovery (% monomer)	Aggregate (% dimer/HMW)	HCP (ppm)	Protein A (ppm)
			Load		3.9/0.54	1260	1.2
1	2	260	3	63.1	2.6/0.09	141	0.009
		260	6	77.5	2.4/0.07	133	0.033
		260	9	82.6	2.4/0.07	134	0.040
2	1	290	3	78.3	3.0/0.13	161	0.007
		290	6	86.3	3.0/0.12	165	0.042
		290	9	88.7	3.1/0.12	170	0.060
3	3	320	3	83.4	4.3/0.20	171	0.003
		320	5	87.7	4.4/0.19	181	0.039
		320	7	90.2	4.5/0.19	185	0.067

Table 3

ANALYSIS DATA OF THE COLLECTED ELUATE PEAKS FROM TOYOPEARL SULFATE-650F

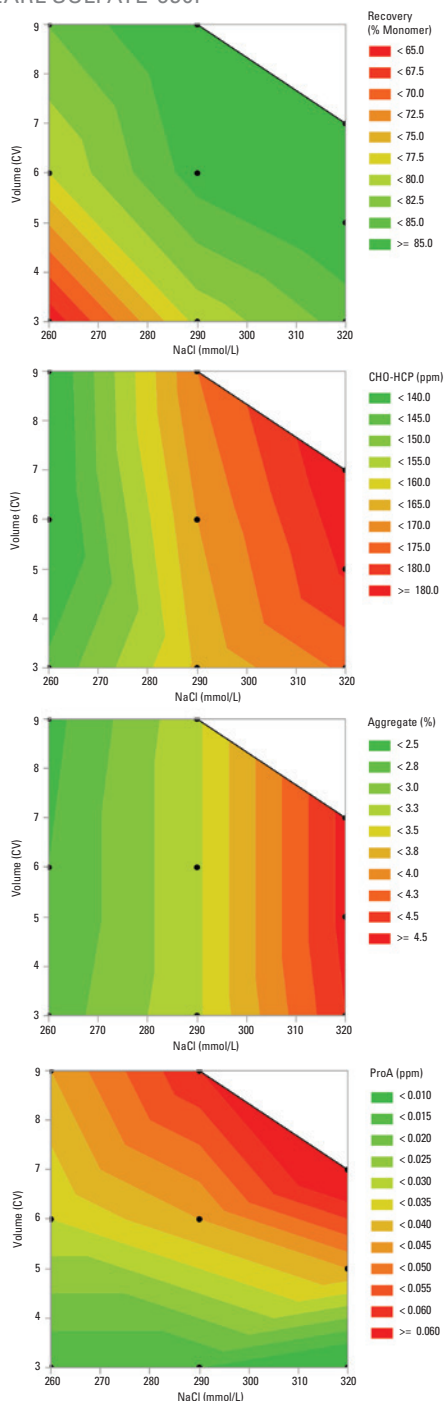



Figure 5

INTEGRATED PEAK DATA FROM SEC COLUMN

Impurity	ProA eluate	Sulfate eluate
Dimer %	3.9	2.4
HMW %	0.54	0.07
HCP (ppm)	1260	134
ProA (ppm)	1.2	0.6040

Table 4

COLLECTED MONOMER PEAKS ANALYZED BY SEC COLUMN

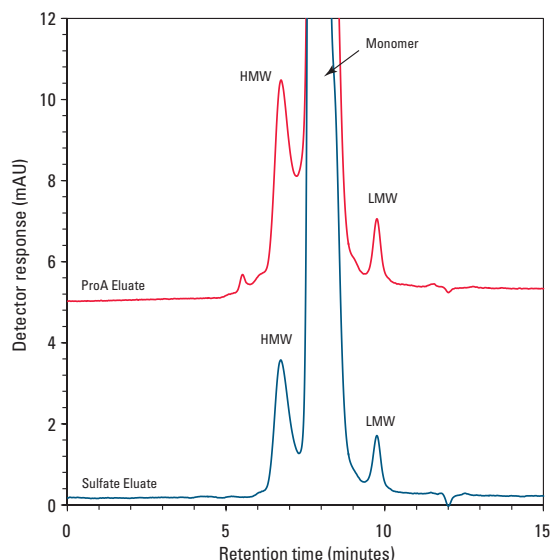

Column: TSKgel G3000SW<sub>XL</sub> (7.8 mm ID x 30 cm)  
 Mobile phase: 0.1 mol/L phosphate, 0.1 mol/L sodium sulfate, 0.1% sodium azide, pH 6.7  
 Gradient: isocratic  
 Flow rate: 1.00 mL/min  
 Detection: UV @ 280 nm (mAU)  
 Temperature: 25° C  
 Injection vol.: 10 µL (100 µg total IgG)  
 Sample: ProA or Sulfate eluate as indicated

Figure 6

Figure 6 shows data from the SEC analysis of the eluate pool at 260 mmol/L NaCl, 9 CV volume. Data shows there is a reduction in aggregate content (in particular HMW impurities), relative to the collected IgG<sub>1</sub> eluate peak material from the protein A resin eluate peak.

The peaks from the SEC column were analyzed for high molecular weight, HCP and protein A ligand content. Table 4 shows that after passing through the TOYOPEARL Sulfate-650F resin, the collected IgG<sub>1</sub> peak has significantly reduced amounts of HMW, HCP and protein A ligand. This suggests that TOYOPEARL Sulfate-650F resin can effectively remove and reduce impurities of IgG<sub>1</sub>.

CONCLUSIONS

The TOYOPEARL Sulfate-650F resin offers a high dynamic binding capacity (>120 mg/mL-resin) with DBC maxima at pH 4.8, 150 mmol/L NaCl and pH 5.2, 100 mmol/L NaCl. With elution at pH 5.2, recovery and impurity removal (aggregate, HCP, leached protein A) is optimal. In fact, analyzed data of the collected IgG<sub>1</sub> monomer peak from the TOYOPEARL Sulfate-650F resin column showed that its purity was significantly improved with an acceptable amount of HMW proteins and HCP while nearly no protein A ligand was detected in the collected IgG<sub>1</sub> peak. By selecting this strong cation exchange resin as a step after mAb post-protein A purification, only a minimal adjustment to pH or salt concentration to the sample is needed.