





# LENS<sup>3</sup>

Multi-Angle Laser Light Scattering

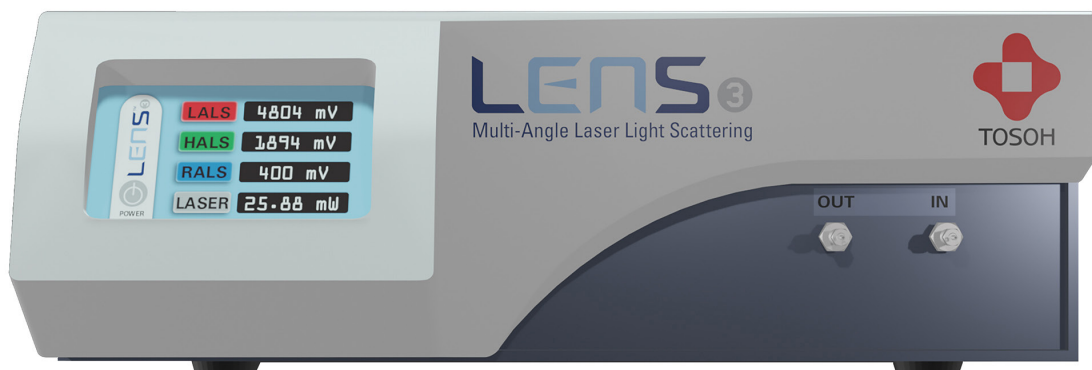

## A Revolutionary Technology for Macromolecular Characterization

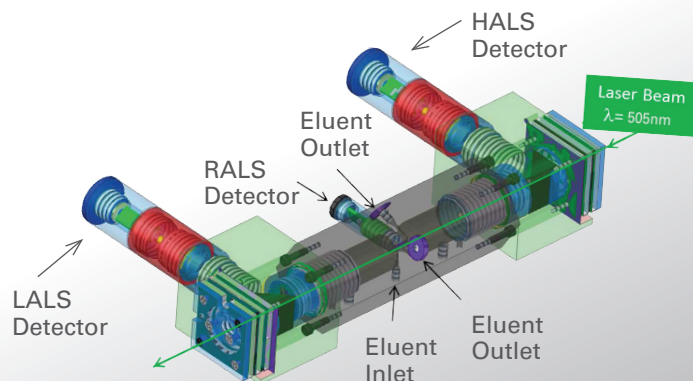

### HIGHLIGHTS

#### High Sensitivity

- Unique patent pending optics design
- Green light source: ~2.7x higher scattering intensity

#### Extended $R_g$ Measurements

- A new calculation method using the Angular Dissymmetry Plot
- Patent pending revolutionary normalization procedure
- Better signal-to-noise ratio with improved electronics and optics

#### SEC/HPLC/UHPLC Compatibility

- High acquisition rate (up to 50 Hz)
- Allows usage of semi-micro SEC/HPLC columns and narrow-bore UHPLC columns

#### Powerful and Intuitive Software

- Simultaneous multi-method execution and analysis
- Simple calibration procedures
- MW and  $R_g$  in a few clicks
- dn/dc and UV extinction coefficient measurement

## Specifications

Number of measurement angles	3
Position of the measurement angles	LALS (10°) RALS (90°) HALS (170°)
Cell geometry	proprietary conical flow path (single inlet, dual outlets)
Laser source type	diode
Laser power	30 mW
Laser wavelength	505 nm
Laser temperature control	yes
Wetted material	teflon, PEEK, glass, stainless steel
Maximum flow rate	2 mL/min
Inlet position	front or side
Baseline noise (in THF @ 1 mL/min)	< 1 mV
Baseline drift (in THF @ 1 mL/min)	< 1 mV / 30 min
MW range	< 200 to 10 <sup>7</sup> Da
R <sub>g</sub> range	< 2 nm to > 50 nm (in progress)
Acquisition rate	< 1 to 50 Hz, user selectable
A/D board channels / resolution	8 channels / 24 bits
Dynamic range	+/- 10 V
Analog inputs	RI, UV and start signal
Connection to PC	ethernet
Dimensions	36.5 (W) x 48.5 (D) x 13 (H) cm = 1.4" x 1.9" x 0.5"
Weight	16 kg = 35 lbs
Intellectual property	PCT/US19/12090: Light Scattering Detectors and Sample Cells for the Same PCT/US19/12095: Light Scattering Detectors and Methods for the Same

### Features and Benefits:

#### MW determination options:

- Direct measurement using 10°
- Direct measurement using 90°
- No longer necessary are Zimm plot extrapolations using multi-angle measurement

#### R<sub>g</sub> determination options:

- A patent pending method using a novel Angular Dissymmetry Plot (no concentration information necessary)
- Replaces the historic assumption of total isotropic scattering for lower molecular size

#### SECview software Features:

- Simultaneous multi-method execution and analysis
- Multi-point  $dn/dc$  and UV extinction coefficient determination
- Automatic peak detection for conventional column calibration method
- Multi-peak selection and "independent" data processing
- Adaptable multiple-injection overlay platform
- Advanced peak band broadening and inter-detector volume correction algorithms
- Direct access to the raw data signals while offering powerful de-spiking and smoothing options
- Access to data point cursors on chromatograms and derived graphs
- Easy export of raw data to ASCII files and graph/chromatogram to picture file

#### Contact us

Dr. Subin Damodaran  
Product Specialist  
subin.damodaran@tosoh.com  
Tel: +49 1628590978

Lisa Jordan  
Technical Inside Sales Specialist  
lisa.jordan@tosoh.com  
Tel: +49 172 9354607